

*...for people in buildings*



# Evacuating in a lying position

... in the event of  
an emergency

**ESCAPE**

MOBILITY COMPANY

**EVACUATION**<sup>®</sup>



# ESCAPE

## MOBILITY COMPANY

® ...for people in buildings

### Escape-Mattress® SOLO

## Evacuating in a lying position ... in the event of an emergency

### Evacuation aids

If use of the lift is not permitted, then evacuation must take place down the stairs.

With an Escape-Mattress®, a vulnerable person, who needs considerable assistance, can be evacuated in a lying position, safely down the stairs.

We advise having an Escape-Mattress®, readily available for those vulnerable persons present, who cannot be evacuated in a sitting position or when the escape route is unfit for the use of an Escape-Chair®.



### A-label quality

The products which we offer have been tested thoroughly to ensure that they meet our high quality standards. Escape Mobility Company carries an 'A-label quality' logo on all of its own products.



**No two buildings are the same and no two situations are the same. Our expert product specialists will be glad to give a demonstration and provide suitable advice. We attach great value to your choice of the right evacuation aids and therefore offer you this service free of charge!**

*Interim modifications in technical specifications, printing and composing errors reserved.*



The Escape-Mattress® Solo is fastened to the wall. By opening the zip at the underside, the Escape-Mattress® Solo slides out of its cover and is ready for use. Suitable for places where there is limited space such as operating theatres.

### Specifications:

- Supplied with a storage cover
- Can be suspended from the wall
- 4 handles on each side;  
1 long and 1 short towing strap at the head end as well as the foot end
- Protection of evacuee/patient:
  - 2 body flaps - 1 per side
  - Foot bag
  - Headrest



Article Number: 11002

## Escape-Mattress® ACCESSORIES

### Escape-Locker

#### Specifications:

- Dimensions: ca. 130 x 80 x 30 cm (supplied without contents)
- Developed especially for storing the Escape-Mattress® Fold or an Escape-Chair®
- Available in white and green
- Can be fastened to the wall



Article Number: 14002

### Escape-Mattress® Storage bag

#### Specifications:

- Dimensions: ca. 125 x 70 x 22 cm
- This storage bag protects the Escape-Mattress® Fold against dust and dirt.
- Two handles at the sides make it easy to transport this storage bag including the Escape-Mattress® Fold.



**NEW!**

Article Number: 11010

# Escape-Mattress® FOLD

# Escape-Mattress® COMPACT

The Escape-Mattress® Fold can be folded double so that it is easy to move and store. Once unfolded, this evacuation mattress is immediately ready for use. This mattress can be stored, hygienically and dust-free, by using one of the accessories.

**Optional accessories:**

- Escape-Locker
- Escape-Mattress® Storage bag

**Specifications:**

- Can be stored folded double
- 4 handles on each side; 1 long and 1 short towing strap at the head end as well as the foot end
- Protection of evacuee/patient:
  - 4 body flaps - 2 per side
  - Foot bag
  - Headrest



Article Number: 1100010

The Escape-Mattress® Compact is a complete Escape-Mattress® but its smaller dimensions allow it to be used flexibly. The Escape-Mattress® Compact is rolled up and stored in the storage bag included. This Escape-Mattress® Compact is more suitable, so that it can be taken to different situations where an evacuation aid is required.

**Specifications:**

- Can be rolled up, comes in a storage bag
- 4 handles on each side; 1 long and 1 short towing strap at the head end as well as the foot end
- Protection of evacuee/patient:
  - 2 body flaps - 1 per side
  - Foot bag
  - Headrest



Article Number: 11000

## Escape-Mattress® COMPARISON TABLE

	Escape-Mattress® Fold	Escape-Mattress® Solo	Escape-Mattress® Compact	Escape-Mattress® Healthcare	Evacuationsheet
Free demonstration and expert advice in advance	✓	✓	✓	✓	✓
Ready for use within a few seconds	✓	✓	✓	✓	✓
Direct evacuation from the bed	not applicable	not applicable	not applicable	✓	✓
Dimensions (length x width)	ca. 193 x 68 cm	ca. 193 x 75 cm	ca. 188 x 68 cm	ca. 198 x 88 cm	ca. 200 x 90 cm
Mattress thickness	ca. 8 cm	ca. 8 cm	ca. 6 cm	ca. 15 cm	not applicable
Including storage bag or case	Optional	✓	✓	-	-
Fits in an Escape-Locker	✓	-	-	not applicable	not applicable
Maximum carrying capacity	150 kg	150 kg	150 kg	150 kg	150 kg
Space-saving	Folding double	Wall mounted	Can be rolled up	not applicable	not applicable
Weight	ca. 7.2 kg	ca. 9.5 kg incl. cover	ca. 7.5 kg incl. bag	ca. 18 kg	ca. 650 gram
Adjustable headrest with immobilisation band	✓	✓	✓	not applicable	not applicable
Warranty	2 years	2 years	2 years	2 years	2 years



## Escape-Mattress® STRETCHER

This Escape-Mattress® Stretcher replaces the standard mattress on an ambulance stretcher. If it is not possible to get to the patient with a stretcher (e.g. because of stairs or other obstacles) this Escape-Mattress® Stretcher can be folded out into a mattress which can be dragged by towing straps at the ends and body flaps for immobilising the patient. After that the mattress with the patient can be lifted onto the stretcher.

### Specifications:

- Dimensions: ca. 180 x 47 cm
- Mattress thickness: ca. 8 cm
- Maximum carrying capacity: 150 kg
- 4 handles on each side;
- 2 adjustable towing straps
- Headrest for fastening onto the stretcher or directly onto the Escape-Mattress®
- Protection of the patient:
  - 2 body flaps for immobilisation
  - Loops for the safety belt of the stretcher
  - Foot bag



Article Number: 1100100

## Evacuationsheet


The Evacuationsheet is positioned permanently between the mattress and the bed frame. In the event of an emergency, the safety straps of the evacuationsheet are secured over the mattress and the patient. The patient is immobilised safely for transporting down the stairs. The Evacuationsheet is used frequently in hospitals and nursing homes.



Article Number: 1100210

### Specifications:

- Dimensions: ca. 200 x 90 cm
- 100% thermally-fixed polyester
- Washable
- Also available in other sizes

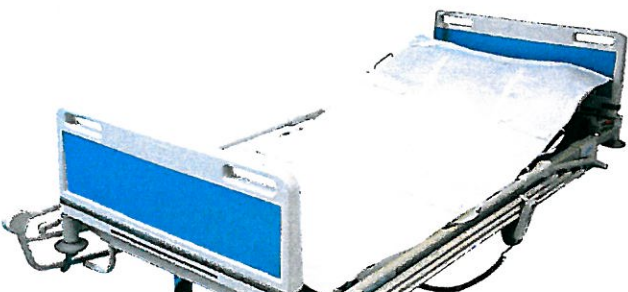

## Escape-Mattress® HEALTHCARE


The Escape-Mattress® Healthcare is an **all-in-one evacuation solution** for use in health care.

This Escape-Mattress® has the appearance and the comfort of an 'ordinary' mattress in a healthcare environment with the difference that this Healthcare evacuation mattress can be deployed immediately, without using other aids, in the event of an evacuation.

In the event of an emergency, you open the zips on the sides of the mattress to pull the body flaps and the towing straps out of the mattress. The body flaps on the sides and at the foot end of the mattress are folded over the bed clothes and the patient. In this way the patient is immobilised and the 'ordinary' mattress is transformed into an evacuation mattress.

### Specifications:

- Standard dimensions: ca. 198 x 88 cm (other mattress sizes are possible)
- Mattress thickness: ca. 15 cm
- Maximum carrying capacity: 150 kg
- 4 handles on each side
- 2 towing straps
- Can be used as an 'ordinary' mattress
- Resistant to bodily fluids
- Anti-bacterial membrane in mattress filling
- Protection for the patient: 3 body flaps for immobilisation

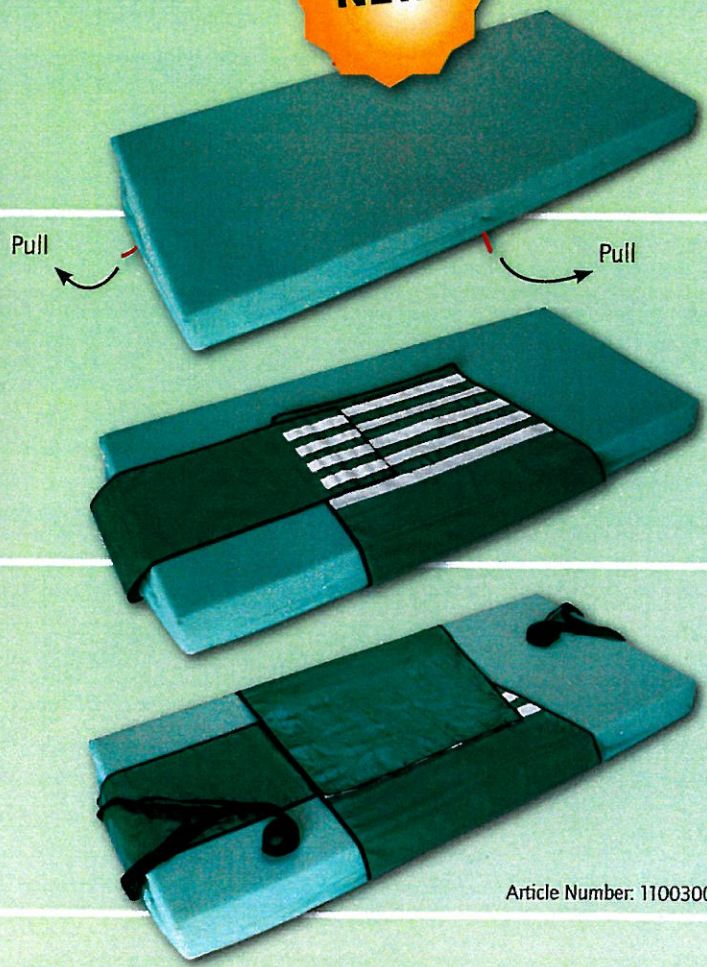

Article Number: 1100300